
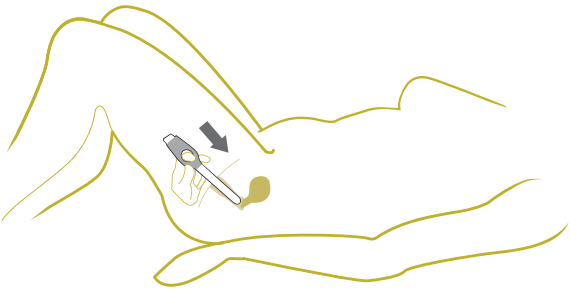


Read these instructions carefully before using the Delphi Screener

The Delphi Screener enables you to collect a fluid sample from your cervix and vagina in a simple and pain free way. It has been designed to make sure that the top will get close to the cervix with insertion in the vagina.


The top contains a small amount (1 teaspoon) of sterile water. By pressing down the button at the back, the sterile water will spread around the cervical area. When the button is released, the water will flow back into the Screener with fluid from the cervix and vagina. After removing the Screener from the vagina the sample is then transferred into a vial. A laboratory will test the sample.

Do not use the Screener when you have your period (i.e. when menstruating) or if you have never had intercourse. If you use vaginal medication or irrigation, you must stop using this 2 days before using the Screener. You can continue to use a vaginal ring, condoms or water-based lubricants. So far, there has been no experience of use with pregnant women, however due to the soft and rounded shape of the Screener, no health or safety risks are expected.

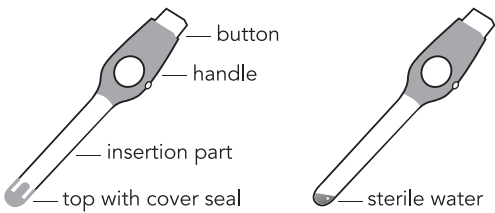

Preparation

- 1 Wash your hands.
- 2 Take out the vial, remove the cap and keep both nearby.
- 3 Open the pouch with the Screener until it is completely open.
- 4 Take the Screener by holding the blue handle.
- 5 Pull the cover seal off the top of the Screener. Some water drops may be visible.

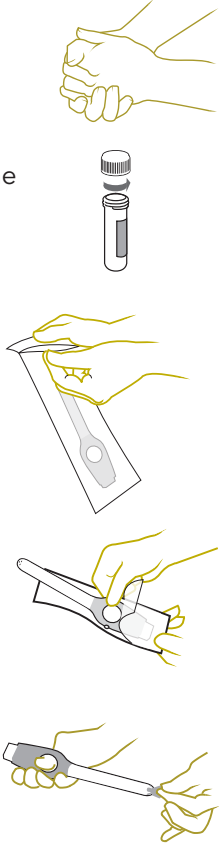

Choose your position


- 6 **A: Lying down**  
Cover a cushion with a towel. Place your buttocks on the cushion. Raise your knees and spread your legs.
- B: Sitting down**  
Sit on a toilet seat with your legs spread.




How to do the self collection

You may have some fluid coming out of your vagina after sampling. This is not a problem. Only a small amount of sample in the vial is needed for the test to be done. You may find some blood or mucus on the Screener - this is normal and will not influence the results of your test.

- 7 Hold the Screener by the blue handle as shown in the picture and use your other hand to spread your labia (vaginal lips). Insert the top of the Screener slowly and gently into the vagina, tilting it backwards, as deep as possible until you feel resistance.
- 8 Hold the Screener whilst pressing the white button fully down.
- 9 Once the button is fully down, count to 3.
- 10 Let go of the button and take the Screener out of your vagina slowly, turning it as you do so. The flushed fluid will then flow back into the Screener. Sometimes the white button does not return immediately to its initial position when taking the Screener out. This is OK.

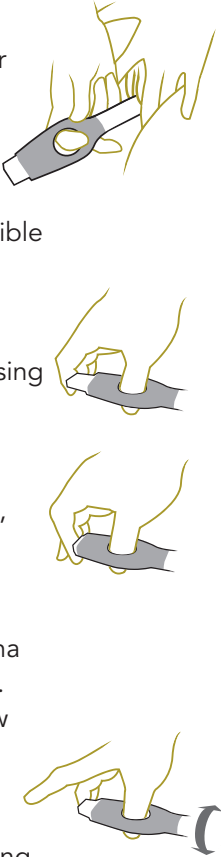

- 11 Insert the Delphi Screener into the vial with the top halfway down. Keep the vial and Screener in an upright position. Press the button of the Screener fully down. The sample will then be released into the vial. The sample will look cloudy and may contain mucus, this is normal.

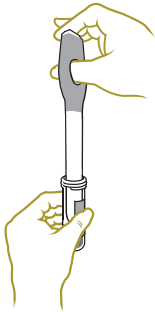


- 12 Take the Screener out of the vial. Place the cap on the vial and make sure it is screwed tightly.




## What to do after taking a sample

13 Check if the vial is correctly labeled with your name.

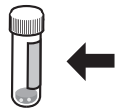


14 Put the Screener and vial with sample in the plastic bag and close both seal strips.

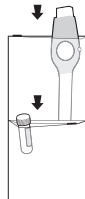

15 Wash your hands again.



16 Hand over the plastic bag to your doctor (-assistant).

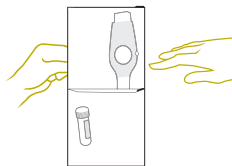

RETURN YOUR SAMPLE WITHIN 24HRS AFTER COLLECTION TO GUARANTEE THE QUALITY OF THE SAMPLE

## General information


You cannot use the Screener after its best before date. This is indicated on the pouch. (🕒 ....)

If the pouch is damaged or opened, the sterility of the Screener cannot be guaranteed and a new one must be used.

Please visit [www.delphiscreener.asia](http://www.delphiscreener.asia) for an animation video on the use of the Delphi Screener, for additional information on the Screener or for User Instructions in Chinese and Malay language.



STERILE R

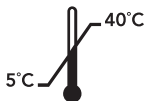

Delphi Bioscience B.V.  
Glashorst 126  
3925 BV Scherpenzeel  
The Netherlands

Date of last revision: February 2012

-210018-

### Disclaimer

Delphi Bioscience cannot accept responsibility in any way for careless or incorrect use of this device. If, despite these instructions, you have any questions or doubts about the use of the device, please contact your doctor.

## User Instructions Delphi Screener

